

bifour
FIAPresto™

ISO 9001 FSC CE IVD

NEW ARRIVAL

Immunofluorescence Analyzer



On site report

Built-in printer

Smart

Auto print & Auto coding

Reliable

*Maintenance free
The selective Lanthanides molecular
was used to label antibody/antigen.
which can reduces the background
noise and improve the sensitivity.*

User friendly interface

7 inch touch screen

Specifications

295×245x100mm
1.7kgs
100-240V AC,50-60Hz

QC/ Calibration

System check card & Controls

Throughput

Up to 20 tests/hour

Interface


7"touch screen
Built-in thermal printer
USB port /Ethernet/SD
card slot
LIS/HIS compatible

Memory

100,000 results

Simple Test Procedure

- STANDARD TEST mode
Automated procedure for single test

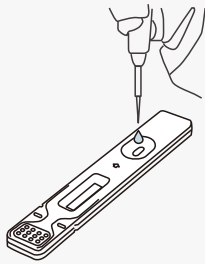

Add the sample into a cassette.




Insert the cassette to start test processing.




Read the result when completed.

- QUICK TEST mode
High-throughput procedure for multiple tests

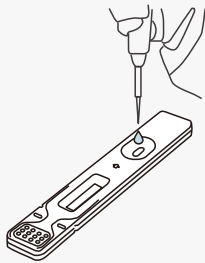

Add the sample into a cassette.

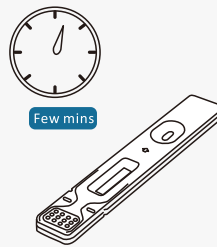

Leave the cassette on the table. Wait for the required time.



Insert the cassette to read the result within 10s.

FIAPresto™ Parameters

Covid-19

- ◆ SARS-CoV-2 Rapid Antigen Test *
- ◆ SARS-CoV-2 & FLU A/B Antigen Combo Test *

Renal Injury Marker

- ◆ NGAL ◆ MAU *
- ◆ β 2-MG *

Bone Metabolism

- ◆ PTH * ◆ CT *
- ◆ 25-OHVD

TDM

- ◆ VM ◆ MTX ◆ VPA
- ◆ CARB ◆ PHT ◆ FK506 *
- ◆ DIG ◆ PTX
- ◆ MPA ◆ THEO

*coming soon

Cardic Marker

- ◆ cTnl ◆ D-Dimer
- ◆ Myo ◆ CK-MB
- ◆ BNP * ◆ LP-PLA2
- ◆ cTnT * ◆ NT-pro BNP
- ◆ HFABP ◆ cTnl/Myo/CK-MB
- ◆ D-Dimer/cTnl/NT-pro BNP *

Inflammation

- ◆ PCT ◆ SAA
- ◆ IL-6
- ◆ PCT&IL-6 *
- ◆ CRP&SAA *
- ◆ hsCRP&CRP

Hepatic Fibrosis

- ◆ CG

Tumor Markers

- ◆ tPSA
- ◆ fPSA
- ◆ AFP *
- ◆ CEA *
- ◆ PGI *
- ◆ PGII *

Thyroid

- ◆ T3
- ◆ T4
- ◆ fT3
- ◆ fT4
- ◆ TSH

Auto-immune

- ◆ IgE

Gonad

- ◆ β -HCG ◆ TES *
- ◆ PROG ◆ E2 *
- ◆ PRL ◆ SHBG *
- ◆ LH
- ◆ FSH
- ◆ AMH *
- ◆ Cortisol *

Anemia

- ◆ Fer

Diabetes

- ◆ CysC
- ◆ HbA1c

Heavy Metal

- ◆ Cd ◆ Pb